



# geoFusion Green Laser

This laser has been Tri-Vortex treated. The special design allows for continuous use with no overheating and this is ideal since no chaos is transferred into the body or liquids.

Use the Green Laser for **PAIN RELIEF** (also for meridian therapy and energy work).

Use the Green Laser to **TREAT WATER** (for better hydration).



## How to Use the geoFusion Green Laser :

- Treat your drinking water or liquids by shining the laser into the liquid for about 1 minute – from top of container onto the surface of the liquid it greatly enhances the quality and taste of the liquid.
- If used for pain relief, shine the laser onto pain area for approx 2 minutes – repeat a few times a day if required.
- Other lasers can NOT do what this one can – if it is not marked "geoFusion", it is NOT recommended for treating of water and for pain relief.

### Tri-Vortex treated –

The focused laser beam is what carries the subtle energy and makes this Green Laser the premium geoFusion product. Many people have reported fast relief from discomfort, soreness, swelling and pain.

### WARNING –

This is a powerful laser and should NOT be pointed into the eyes!



Manufactured by : Tri-Vortex SA,  
PO Box 645, Lonehill, 2062,  
South Africa

# geoFusion Green Laser

## INTRODUCTION:

This is a high quality Diode Pumped Solid State Laser. Use the laser to point at any screens or objects - as far as 500m.

## USER SAFETY PRECAUTIONS :

Under normal operation, it is safe to the human body. Do not point the beam into the face or eyes - it can cause eye injury !

## BATTERY INSTALLATION:

Unscrew the end cap from the casing. Insert 1 x CR123A battery with the + terminal facing towards the cap. Using alkaline batteries is recommended.

## OPERATION :

Press the push button in the cap to turn the laser on. Press the button once more to turn the laser off. Operating temperature : 12-38 degC

## LIMITED WARRANTY :

The manufacturer warrants this product to be free of defects in workmanship and materials for a period of one year after purchase. This warranty (excl battery) is solely limited to the repair and replacement of original parts. All other costs are for the owner. This warranty does not cover any defects, damage and deterioration due to misuse, alteration or negligence. Tampering with any of the labels will void all warranties.


## SPECIFICATIONS:

Wavelength - 532nm  
 Battery - (3V Li) - 1 x CR123A  
 Output Power - < 5mW CLASS IIIA

Approval :



KEEP AWAY FROM CHILDREN

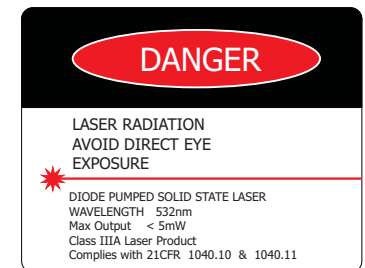

Made in USA